

Self-assembled silane monolayers: A step-by-step high speed recipe for high-quality, low energy surfaces

M. Lessel,[†] O. Bäumchen,^{†,‡} M. Klos,[†] H. Hähl,^{†,¶} R. Fetzer,^{†,§} R. Seemann,[†] and K. Jacobs^{*,†}

Department of Experimental Physics, Saarland University, D-66041 Saarbrücken, Germany

E-mail: k.jacobs@physik.uni-saarland.de

Abstract

Silanization bases on the adsorption, self-assembly and covalent binding of silane molecules onto surfaces, resulting in a densely packed self-assembled monolayer (SAM). Following standard recipes, however, the quality of the monolayer is often variable and therefore unsatisfactory. The process of self-assembly is highly affected by the chemicals involved in the wet-chemical process, by the ambient parameters during preparation such as humidity or temperature and by possibly present contaminants (or impurities). Here, we present a reliable, efficient, and wet-chemical recipe for the preparation of ultra-smooth, highly ordered alkyl-terminated silane SAMs on Si wafers.

Introduction

SAMs prepared on smooth surfaces like Si wafers exhibit extraordinary properties such as chemical homogeneity, ultra low surface roughness and

controlled wettability.^{1–3} The latter can be varied in wide ranges, depending on the end group of the silane. Here, we focus on alkyl-terminated silanes resulting in surfaces of very low surface energy. Silane layers in particular are mechanically robust⁴ and thermally stable up to at least 250 °C.^{5,6} These properties render silane-coated substrates ideal model surfaces to study a wide range of physical, chemical and biological phenomena such as adhesion,^{7–9} adsorption,^{10–13} friction^{14–17} or nanofluidics of thin liquid films.^{18–23} Alongside fundamental research, also technical applications benefit from their unique properties: SAMs act *e.g.* as lubrication layers in micro-electro-mechanical-systems (MEMS)⁵ and as coatings in microfluidic devices.^{24,25}

The recipe presented here focuses on the widely-used CH₃-terminated alkylsilanes with three different chain lengths determined by the number of carbon atoms (12, 16 and 18) in the backbone: dodecyl-trichlorosilane (DTS, C₁₂H₂₅Cl₃Si), hexadecyl-trichlorosilane (HTS, C₁₆H₃₃Cl₃Si) and octadecyl-trichlorosilane (OTS, C₁₈H₃₇Cl₃Si), illustrated in Figure 1.

The formation of self-assembled monolayers was first reported by Zisman^{27,28} in the late 1940's. From the 1980's onwards, Sagiv and others^{2,29–31} studied the underlying formation process and characterized surface properties like surface energy, layer thickness and roughness on alkylsilane SAMs. Meanwhile, a variety of experimental techniques have been applied to study

^{*}To whom correspondence should be addressed


[†]Department of Experimental Physics, Saarland University, D-66041 Saarbrücken, Germany

[‡]Present address: Department of Physics and Astronomy, McMaster University, Hamilton, ON, Canada L8S 4M1


[¶]Present address: University of Zurich, Institute of Physical Chemistry, Winterthurerstrasse 190, 8057 Zurich, Switzerland

[§]Present address: Karlsruhe Institute of Technology, Institute for Pulsed Power and Microwave Technology, 76344 Eggenstein-Leopoldshafen, Germany

DTS

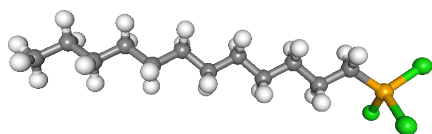

HTS

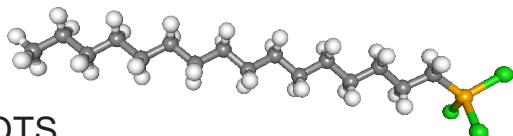

OTS

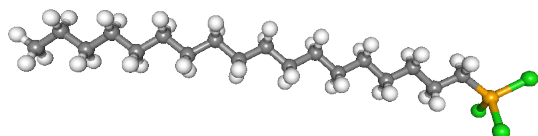


Figure 1: Ball and stick representation of alkylsilane molecules (white: hydrogen, gray: carbon, orange: silicon, green: chlorine). (Models illustrated using BALLView.²⁶)

SAMs: atomic force microscopy,³² ellipsometry,^{29,30} contact angle measurements,^{33,34} low frequency infrared spectroscopy,³⁵ near edge X-ray absorption spectroscopy³³ or X-ray reflectometry.³⁶ Most of these studies followed a wet-chemical preparation procedure, going back to the original work of Zisman.^{27,28} It turns out that the quality of the silane SAM does not only depend on the recipe itself and on the quality of the used chemicals, but also on humidity, temperature and contaminants present during preparation.^{33,34,37,38} In the following, we present a quick and easy-to-follow recipe to prepare silanized Si wafers of reproducible quality, *i. e.* of low roughness, high ordering, high water contact angle and low contact angle hysteresis.

Preparation

Alkylsilanes are offered from several suppliers in very high purity. For our study, DTS (> 99 %), OTS (> 90 %) and HTS (95 %) were used (*c.f.* table 1). Flakes or filaments in the solution indicate aggregations of silane molecules and are typically caused by contact of the solution with water or humid air, inducing cross-linking of the silanes.³⁴ If this is the case, the solution should not be used. Silane solutions are packed under inert gas atmosphere from to supplier to prevent

such cross-linking. We therefore recommend to use a new package of silane for each silanization. A complete list of necessary chemicals is given in Table 1.

The coating process consists of four steps: (i) cleaning of all parts involved in the preparation, such as beakers and tweezers, (ii) cleaning of the Si wafers, (iii) preparation of the solution for the silanization and (iv) the coating procedure itself⁴⁰

Cleaning Procedure

A clean wafer surface is a precondition to achieve a covalent binding of the silane to the natural, amorphous SiO₂ layer of the Si wafer.⁴¹ Typically, organic residues left over from the polishing process of the silicon wafer have to be removed. Here, peroxymonosulfuric acid ('piranha solution'),⁴² consisting of 50 % H₂SO₄ and 50 % H₂O₂ has proven to be suitable as it does not roughen the wafer surface. For details of the chemicals see Table 1.

All steps (i) - (iv) were performed in a dust-free environment, at least a class 100 clean room or equivalent environments (*e. g.* a laminar flow work bench). All beakers⁴³ and tweezers to be used during the silanization process should withstand the cleaning procedure. Containers and tweezers made of Teflon[®] or other PTFE polymers have proven to be suitable.

For step (i), all beakers are filled with piranha solution. Then, the tweezers are given into the acid. To ensure complete cleaning, the beakers (together with the tweezers) are tilted and then slowly rotated by hand (5-10 min), so that the piranha solution can reach the bezel of the beaker. Subsequently, beakers and tweezers are rinsed three times in hot deionized (DI-) water (18.2 MΩcm at 25 °C and total-organic-carbon less than 6 ppb). This routine is repeated twice, each with fresh piranha solution and fresh DI-water. To prevent new contamination, the tweezers are kept in the cleaned containers.

For step (ii), in the first PTFE beaker, one Si wafer (with native SiO₂,⁴⁴ *e.g.* in our case purchased from Si-Mat Silicon Materials, Kaufering, Germany) is immersed (polished side up) into freshly prepared piranha solution straight from the dust-tight package of the wafer manufacturer. The

formation of bubbles on the wafer surface is an indication of the acid's reactivity. A periodical, careful waving of the beaker removes bubbles from the wafer surface and stirs the liquid to ensure uniform oxidation of any organic residues on the surface. After 30 min, the wafer is removed from the piranha solution and placed into the second PTFE beaker filled with hot DI-water. Promptly, the water is replaced by new hot DI-water. The wafer remains for about 10-20 min in the new water. Then, the liquid is replaced every 20 min, at least 3 times. This procedure is necessary to completely remove reaction products left over on the wafer surface. After the piranha treatment, the Si wafer exhibits a very high surface energy and is thus susceptible to contamination. Therefore, the wafer should be kept in the DI-water until needed. While the first wafer is immersed in the DI-water a second wafer can be cleaned with piranha according to step (ii).

Preparation of the Silane Solution

For step (iii), the third PTFE beaker, cleaned in step (ii), is thoroughly dried with an inert gas *e. g.* N₂ (5.0 grade).⁴⁵ The following quantities suffice for immersing one 4 inch wafer at a time into a suitable beaker, requiring about 50 ml of silane solution.

The basic solvent is bicyclohexyl. In 50 ml of bicyclohexyl, 30 drops (approx. 1.5 μ ml) of carbon tetrachloride CCl₄ and 15 drops (approx. 0.75 μ l) of the silane are added. As CCl₄ prevents the silane from forming micelles,³⁴ the solution is lucid. The appearance of white flakes or filaments indicates an uncontrolled cross-linking of silane molecules as described above. In this case, the solution should not be used. We achieved the presented results at relative humidities below 45 % at typical room temperatures (around 25 °C) and with freshly prepared silane solutions.

Coating Procedure

Despite the negative implications of water for the silane solution, the presence of water is also necessary for a certain step of the SAM preparation: Silberzahn et al.,⁴⁶ Angst and Simmons,⁴⁷ Tripp et al.⁴⁸ and Wang et al.³⁷ highlighted that a thin layer (3 nm) of water⁴⁹ adsorbed on the Si wafer

surface improves the quality of the resulting SAM drastically. Implementing this aspect in step (iv), the cleaned Si wafer is taken out of the hot DI-water bath and then dried in a N₂ jet. When no water is visible anymore on the surface, it is exposed shortly (half a second—in order to prevent the formation of droplets) to the vapor of a boiling DI-water bath until a breath figure appears, which should evaporate again within a second. Immediately after the breath figure vanishes, the wafer is put for 15 min into the silane solution, which should completely cover the polished side of the wafer. Afterwards, the wafer is taken out of the silane solution and flushed with a jet of chloroform out of a wash bottle to remove unbound silane molecules.

Especially on humid days, a white 'haze' can build up on the silane solution for the same reason mentioned before: silane molecules cross-link in the presence of water.³⁴ If this is the case, the surface of the solution should be skimmed with an extra piece of a cleaned Si wafer. A similar haze can be observed sometimes on the silanized Si wafer. It can be easily removed by thoroughly rinsing the wafer with chloroform.

Before placing the wafer back into the silane solution, all chloroform residues have to be removed by drying the entire wafer (both sides) with a jet of clean, dry N₂. Repeating the process of rinsing with chloroform and bathing in silane solution enhances the quality of the SAM. Typically, three iterations are sufficient. With some experience, the quality of the SAM can already be judged by the way the chloroform and the silane solution bead off the surface: If there is a smooth, straight three phase contact line between liquid, solid and air, with a high contact angle and a low hysteresis (which can be judged by the trained eye), the silanization was successful and a homogeneous SAM can be expected.

In the following section, the resulting SAM surfaces are characterized in terms of their wetting properties, layer thickness and the surface roughness on different lateral length scales. A structural analysis using X-ray reflectometry of the SAM resulting from this procedure is presented in Ref.⁵⁰

Characterization

Water contact angles and surface energy

Following the preparation procedure presented in this study, advancing water contact angles of $114(2)^\circ$ and a contact angle hysteresis below 10° can reliably be achieved for DTS, HTS and OTS (*c.f.* Table 2). The contact angles were obtained with fresh DI-water using the sessile drop method (OCA 20, Dataphysics Instruments GmbH, Filderstadt, Germany). The surface energy of the silane SAM can be approximated using the Good Girifalco relation,⁵¹ if the contact angles of two apolar liquids are known. Using 1-bromonaphthalene ($> 95\%$, Fluka) and bicyclohexyl (99% , Acros Organics), contact angles on OTS of $40(4)^\circ$ and $60(4)^\circ$, respectively, were achieved. The results for the surface energy of the different silane SAMs are listed in Table 2.

Structural analysis

For the structural analysis of the SAMs, X-Ray reflectometry (XRR)⁵⁰ was applied. By XRR, a layer thickness of $2.17(1)\text{ nm}$ ⁵⁰ was measured for the OTS tail-group, which is in excellent agreement with the theoretically calculated length of the OTS tail in *all-trans* configuration of 2.18 nm ,⁵⁰ indicating an upright orientation. The results for DTS and HTS are compiled in Table 2 and are detailed in Ref.⁵⁰

Surface topography and atomic force microscopy (AFM)

The quality of the SAM in terms of surface coverage and roughness is characterized by atomic force microscopy (AFM, a Multimode and a FastScan Icon, both Bruker Nano Inc., Santa Barbara, CA, USA) in TappingmodeTM.

On a $1\text{ }\mu\text{m}^2$ area, the root mean square (rms) roughness of the OTS layer is determined to $0.17(3)\text{ nm}$, that of a bare Si wafer to $0.14(2)\text{ nm}$ (*c.f.* Figure 2 a and b, 515×512 pixel). Figure 3 exemplifies a $(6\text{ }\mu\text{m})^2$ AFM scan (616×616 pixel) of a typical OTS surface: The surface is homogeneously covered by the SAM and - at variance with other publications^{15,52} - no holes or grain boundaries can be detected. Operating the AFM in

TappingmodeTM enables us to record an additional information, the phase lag between the free end of the cantilever and the driving piezo element, which can be taken as a measure for the energy dissipation.⁵³ On a flat surface ($\text{rms} < 1\text{ nm}$), differences of the phase lag typically hint to a chemically heterogeneous surface.⁵⁴ Figure 2 c and d depict the phase signal for the OTS-SAM and a bare Si wafer. On both surfaces, a very low phase lag (below 2°) is recorded, which is an additional indication for the chemical homogeneity of the surface.

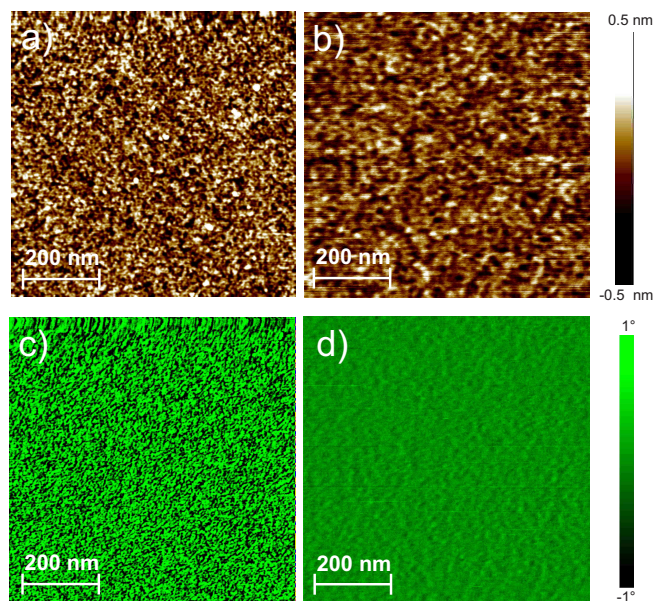


Figure 2: AFM images: Topography of an OTS-SAM (a) in comparison with a bare Si wafer (b) prior to the preparation of the SAM and corresponding phase signal (c) of an OTS-SAM and (d) Si wafer.

A similar high quality surface can be found for HTS and DTS (*c.f.* Table 2 and Ref.⁵⁰). Though the preparation temperature of 25°C is above the ideal coating temperature reported by Brzoska et al.³⁴ for the three types of silanes studied, the rms roughness of the different SAMs are similar (*c.f.* Table 2). The tilt angles of the hydrocarbon chains, however, are different as is detailed in the XRR study by Gerth et al.⁵⁰ (*c.f.* Table 2).

The following section compares the recipe presented here with the most common procedures reported in literature concerning time efficiency of the preparation procedure.

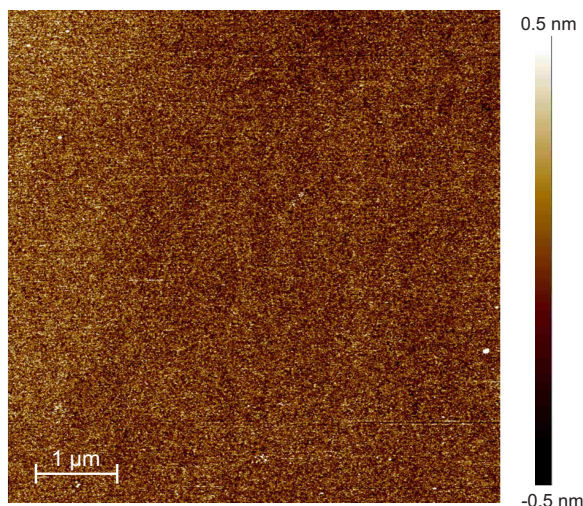


Figure 3: AFM height data of a representative $(6\text{ }\mu\text{m})^2$ scan area (616×616 pixels) obtained from OTS. The image visualizes the large-scale homogeneity and full surface coverage of the SAM.

Comparison

Since the pioneering work of Zisman in 1946,^{27,28} a variety of recipes to prepare monolayers of silanes via self-assembly on SiO_2 substrates were reported.

However, chemical cleaning of the Si wafers always represents the first step of the SAM preparation. Based on either alkaline or acidic solutions, various procedures are reported in literature: Wang et al.³⁷ applied RCA1 and RCA2,⁵⁵ two procedures well established in the semiconductor industry. Sagiv²⁹ used nitric acid, Wasserman et al.³⁰ and Mezger et al.³⁶ as well as this recipe used piranha solution to clean the SiO_2 surfaces. A ozone plasma treatment is reported by Brzoska et al.³⁴ Most recipes involve rinsing of the substrates with DI-water to remove aggressive cleaning chemicals.

An important feature to our experience is the cleaning of the containers and tweezers, which, however, is described only by a few articles, *e. g.* Wang et al.³⁷ and this study.

In terms of the chemical composition of the silane solution, different options can be found in the literature. The most common solvents for the silane solution are: hexadecane,^{30,34} Isopar G³⁷ or - as in this recipe - bicyclohexyl.^{30,32} To avoid the agglomeration of micelles caused by water in the solution, small amounts of either chloroform

or carbon tetrachloride are added.

The time period the Si wafer is exposed to the silane solution differs significantly from recipe to recipe: typical exposure times range from several minutes³⁴ to 10 days.³⁷ Afterwards, the SAMs are rinsed with a solvent (chloroform (CHCl_3)³⁰ or dichloromethane (CH_2Cl_2)^{30,32}) to remove unbound silanes. As described above, the silane coated substrates are rinsed with chloroform several times during the coating process, which improves the quality of the resulting SAM.

Several articles describe different kinds of post-treatment of the prepared SAMs to enhance the cross-linking of the silanes such as heating in an oven³² or storage in a desiccator.³⁷ The recipe presented here features only a rinse in chloroform as the last step. The coated samples can then be stored in a dust-controlled environment in air. As compared to Si wafers, the silanized wafers exhibit a very low surface energy. Physisorbed particles, generated for instance by cutting the wafers, can usually be flushed away easily, *e. g.* by sonicating the samples in polar and apolar ultra-clean solvents like ethanol, acetone or toluene⁵⁶ for a few minutes.

Conclusion

In conclusion, we present a reliable and easy-to-follow recipe that leads to silane SAMs of high quality in terms of a low contact angle hysteresis and a very low surface roughness (comparable to a bare Si wafer), both resulting from the excellent topographical and chemical homogeneity of the SAMs. The characterization of the SAMs includes AFM, contact angle measurements, ellipsometry and X-ray reflectometry.⁵⁰ In addition, it is optimized concerning preparation time and allows the silanization of a complete Si wafer in about 3 h. The throughput can be enhanced, if the cleaning and coating process is simultaneously performed for several wafers. With some practice and sufficient amount of clean beakers, 3-4 wafers can be prepared in one run, which corresponds to an average preparation time of 45 to 60 min per wafer. The recipe has been tested in several labs with different environmental conditions with similar experience and similar quality. With this study,

we like to set a lab standard for commonly used silane coatings on Si wafers at normal lab conditions. The intention is to facilitate a comparison of silane coatings of different studies and different labs. Standard quality checks of a silane SAM should comprise at least an advancing and receding water contact angle measurement as well as the determination of the sample topography of a 1-10 μm -sized area and roughness before and after coating.

Notes and References

- (1) Ulman, A. *Chem. Rev.* **1996**, 96, 1533–1554.
- (2) Schreiber, F. *Prog. Surf. Sci.* **2000**, 65, 151–257.
- (3) Ravoo, B. J.; Reinhoudt, D. N.; Onclin, S. *Angew. Chem. Int. Ed.* **2005**, 44, 6286–6304.
- (4) Ding, J. N.; Wong, P. L.; Yang, J. C. *Wear* **2006**, 260, 209–214.
- (5) Srinivasan, U.; Houston, M. R.; Howe, R. T.; Maboudian, R. *J. Microelectromech. Syst.* **1998**, 7, 252–260.
- (6) Helmy, R.; Fadeev, A. Y. *Langmuir* **2002**, 18, 8924–8928.
- (7) Prime, K. L.; Whitesides, G. M. *Science* **1991**, 252, 1164–1167.
- (8) Prime, K. L.; Whitesides, G. M. *J. Am. Chem. Soc.* **1993**, 115, 10714–10721.
- (9) Faucheux, N.; Schweiss, R.; Lützow, K.; Werner, C.; Groth, T. *Biomaterials* **2004**, 25, 2721–2730.
- (10) Cao, T.; Tang, H.; Liang, X.; Wang, A.; Auner, G. W.; Salley, S. O.; Ng, K. Y. S. *Biotechnol. Bioeng.* **2006**, 94, 167–176.
- (11) Tang, H.; Cao, T.; Wang, A.; Liang, X.; Salley, S. O.; McAllister II, J.; Ng, K. *J. Biomed. Mater. Res.* **2007**, 80A, 885–894.
- (12) Hähl, H.; Evers, F.; Grandthyll, S.; Paulus, M.; Sternemann, C.; Loskill, P.; Lessel, M.; Hüsecken, A. K.; Brenner, T.; Tolan, M.; Jacobs, K. *Langmuir* **2012**, 28, 7747–7756.
- (13) Loskill, P.; Hähl, H.; Kreis, C. T.; Thewes, N.; Bischoff, M.; Herrmann, M.; Jacobs, K. *Langmuir* **2012**, 28, 7242–7248.
- (14) Tsukruk, V.; Bliznyuk, V. N. *Langmuir* **1998**, 14, 446–455.
- (15) Flater, E. E.; Ashurst, W. R.; Carpick, R. W. *Langmuir* **2007**, 23, 9242–9252.
- (16) Chandross, M.; Lorenz, C. D.; Stevens, M. J.; Grest, G. S. *Langmuir* **2008**, 24, 1240–1246.
- (17) Booth, B. D.; Vilt, S. G.; Lewis, J. B.; Rivera, J. L.; Buehler, E. A.; McCabe, C.; Jennings, G. K. *Langmuir* **2011**, 27, 5909–5917.
- (18) Pit, R.; Hervet, H.; Léger, L. *Phys. Rev. Lett.* **2000**, 85, 980–983.
- (19) Cottin-Bizonne, C.; Jurin, S.; Baudry, J.; Crassous, J.; Restagno, F.; Charlaix, E. *Eur. Phys. J. E.* **2002**, 9, 47–53.
- (20) Fetzer, R.; Jacobs, K.; Münch, A.; Wagner, B.; Witelski, T. P. *Phys. Rev. Lett.* **2005**, 95, 127801.
- (21) Fetzer, R.; Rauscher, M.; Münch, A.; Wagner, B. A.; Jacobs, K. *Europhys. Lett.* **2006**, 75, 638–644.
- (22) Zhang, Z.; Wang, Z.; Xing, R.; Han, Y. *Surf. Sci.* **2003**, 539, 129–136.
- (23) Zhu, X.; Wang, H.; Liao, Q.; Ding, Y. D.; Gu, Y. B. *Exp. Therm. Fluid. Sci.* **2009**, 33, 947–954.
- (24) Tretheway, D. C.; Meinhart, C. *Phys. Fluids* **2002**, 14, L9–L12.
- (25) Tretheway, D. C.; Meinhart, C. *Phys. Fluids* **2004**, 16, 1509–1515.
- (26) Hildebrandt, A.; Dehof, A. K.; Rurainski, A.; Bertsch, A.; Schumann, M.; Toussaint, N.; Moll, A.; Stockel, D.; Nickel, S.; Mueller, S.; Lenhof, H. P.; Kohlbacher, O. *BMC Bioinformatics* **2010**, 11, 531.

- (27) Bigelow, W. C.; Pickett, D. L.; Zisman, W. A. *J. Colloid. Sci.* **1946**, *1*, 513–538.
- (28) Bigelow, W. C.; Glass, E.; Zisman, W. A. *J. Colloid. Sci.* **1947**, *2*, 563–591.
- (29) Sagiv, J. *J. Am. Chem. Soc.* **1980**, *102*, 92–98.
- (30) Wasserman, S. R.; Tao, Y. T.; Whitesides, G. M. *Langmuir* **1989**, *5*, 1074–1087.
- (31) Banga, R.; Yarwood, J.; A. M. Morgan, B. E.; Kells, J. *Langmuir* **1995**, *11*, 4393–4399.
- (32) Bierbaum, K.; Grunze, M.; Baski, A. A.; Chi, L. F.; Schrepp, W.; Fuchs, H. *Langmuir* **1995**, *11*, 2143–2150.
- (33) Peters, R.; Nealey, P. F.; Crain, J. N.; Himpfel, F. J. *Langmuir* **2002**, *18*, 1250–1256.
- (34) Brzoska, J.; Azouz, I. B.; Rondelez, F. *Langmuir* **1994**, *10*, 4367–4373.
- (35) Tripp, C. P.; Hair, M. L. *Langmuir* **1995**, *11*, 1215–1219.
- (36) Mezger, M.; Schöder, S.; Reichert, R. H.; Schröder, H.; Okasinski, J.; Honkimäki, V.; Ralston, J.; Bilgram, J.; Roth, R.; Dosch, H. *J. Chem. Phys.* **2008**, *128*, 244705.
- (37) Y. Wang, M. L. *Langmuir* **2003**, *19*, 1159–1167.
- (38) Wang, M.; Liechti, K. M.; Wang, Q.; White, J. M. *Wear* **2005**, *21*, 1848–1857.
- (39) It is essential to use non-stabilized H₂O₂.
- (40) For handling the chemicals, local safety rules apply and the use of lab coats, goggles and suitable gloves is mandatory.
- (41) Loskill, P.; Hähl, H.; Faidt, T.; Grandthyll, S.; Müller, F.; Jacobs, K. *Adv. Cool. Interf. Sci.* **2012**, *107*, 179182.
- (42) A compact guide for handling 'piranha solution' can be found *e.g.* at <http://www.nanofab.ubc.ca/process-piranha-etch>.
- (43) For the silanization of a 4 inch wafer, three PTFE beakers of similar diameter,⁵⁷ one 2 l Erlenmeyer flask for the hot DI-water, one graduated cylinder (50 ml) and pipettes for the chemicals of the silanization solution.
- (44) Crystal orientation and type of doping of the Si-wafers do not impact on the silanization.
- (45) Including an oil free blow gun.
- (46) Silberzahn, P.; Leger, L.; Auserre, D.; Benattar, J. J. *Langmuir* **1991**, *7*, 1647–1651.
- (47) Angst, D. L.; Simmons, G. W. *Langmuir* **1991**, *7*, 2236–2242.
- (48) Tripp, C. P.; Hair, M. L. *Langmuir* **1992**, *8*, 1120–1126.
- (49) Asay, D. B.; Kim, S. H. *J. Phys. Chem. B* **2005**, *109*, 16760–16763.
- (50) Gerth, S.; Lessel, M.; Bäumchen, O.; Klimczak, M.; Steinrück, H. G.; Jacobs, K.; Magerl, A. *Langmuir* **2012**, *in Preparation*.
- (51) Good, R. J.; Girifalco, L. A. *J. Phys. Chem.* **1960**, *64*, 561–565.
- (52) Carraro, C.; Yauw, O. W.; Sung, M. M.; Maboudian, R. *J. Phys. Chem. B* **1998**, *102*, 4441–4445.
- (53) Tamayo, J.; Garcia, R. *Appl. Phys. Lett.* **1998**, *73*, 2926–2928.
- (54) Knoll, A.; Krausch, G. *Macromolecules* **2001**, *34*, 4159–4165.
- (55) Kern, W.; Puotinen, D. *RCA Rev.* **1970**, *31*, 187.
- (56) Especially for toluene, we recommend using semiconductor industry-grade quality for its particle-controlled content, for instance LiChroSolv® (Sigma-Aldrich) or Optima® (Fisher Scientific).
- (57) It is convenient to choose the diameter of the beaker only slightly larger than the diameter of the wafer as this allows a complete coverage of the wafer with a minimum amount of piranha solution and silane solution required.

Acknowledgement

We gratefully acknowledge financial support from German Science Foundation DFG under Grant No. Ja 905/3 within the priority program 1164, Grant No. Se 1118/2 and the graduate school GRK 1276 for financial support.

Table 1: Chemicals necessary for silanization.

Chemicals	purity /or special product	CAS-no.	supplier
<u>cleaning:</u>			
hydrogen peroxide	30 %, non-stabilized ³⁹	7722-84-1	Merck
sulphuric acid	96 %, Selectipur [®]	7664-93-9	BASF
<u>solvents:</u>			
bicyclohexyl	99 %	92-51-3	Acros Organics
carbon tetrachloride		56-23-5	local supplier used
chloroform	> 99.8 %, LiChrosolv [®]	67-66-3	Sigma-Aldrich
<u>silanes:</u>			
octadecyl-trichlorosilane	> 90 %	112-04-9	Sigma-Aldrich
hexadecyl-trichlorosilane	95 %	5894-60-0	ABCR
dodecyl-trichlorosilane	puriss. > 99 %	4484-72-4	Fluka

Table 2: Surface properties of silane SAMs.

silane	$H_2O\theta_{adv}$	$\Delta\theta$	surface energy / mNm^{-1}	thickness silane tail (XRR) ⁵⁰ / nm	tilt angle of the silane tail ⁵⁰	roughness (rms) $1\mu m^2$ AFM scan / nm
DTS	114(2)°	5°	26(1)	1.31(3)	22(1)°	0.13(3)
HTS	114(2)°	8°	25(1)	1.87(3)	14(1)°	0.16(2)
OTS	111(2)°	5°	24(1)	2.17(1)	5(1)°	0.17(2)

Table 3: Comparison of the most common silanization recipes.

	container cleaning	substrate cleaning	silanisation solution	finishing treatment
Sagiv ²⁹ (1980)	not reported	30 min (HNO ₃ , NaOH),	no time reported (HCl ₃ , CCl ₄ , alkane)	not reported
Wasserman ³⁰ (1989)	not reported	≈ 1 h (piranha solution)	1-48 h (hexadecane or bicyclohexyl, CCl ₄)	CH ₂ Cl ₂ , CHCl ₃ , ethanol
Brzoska ³⁴ (1994)	not reported	several hours (piranha solution, O ₃ plasma)	between 10 min and 1 h (CCl ₄ , alkane)	not reported
Wang ³⁷ (2002)	48 h	≈ 4 h (CHCl ₃ acetone, ethanol, RCA1, RCA2)	18 h to 19 d (isoprop G, acetone, CHCl ₃)	CH ₂ Cl ₂ , acetone, alcohol, 12 h in a desiccator
Mezger ³⁶ (2008)	not reported	≈ 1 h (UV in O ₂ atmosphere and 10 min in piranha solution)	3 h (hexadecane, CHCl ₃)	hexadecane, toluene, 1 h at 100 ° C
this study	≈30 min (piranha solution)	around 1.5 h (piranha solution)	3x15 min (CCl ₄ , bicyclohexyl), with CHCl ₃ in between	rinsed with CHCl ₃