

Analyzing Single Crystal Time-of-Flight Neutron Data – NOBUGS2002/030

P.F. Peterson,¹ D.J. Mikkelsen,² R.L. Mikkelsen,² A.J. Schultz,¹ J.P. Hammond,¹ A. Chatterjee,¹ and T.G. Worlton¹¹Intense Pulsed Neutron Source, Building 360, 9700 South Cass Ave, Argonne, IL 60439-4814, USA²Department of Mathematics, Statistics and Computer Science,
University of Wisconsin-Stout, Menomonie, WI 54751, USA

(Dated: Oct 19, 2002)

Software for time-of-flight single crystal diffractometer (SCD) data visualization and analysis was originally written to be run in batch mode on VMS-VAX systems. Modern computers and software tools allow a new level of visualization and user interaction. The Integrated Spectral Analysis Workbench (ISAW) is being extended and customized for SCD measurements. To this end, new viewers and operators have been added to ISAW. One of the new operators allows running a program written in another language, such as C or FORTRAN, to integrate tested software into the object-oriented package. This provides a method for preserving earlier software development while making use of modern tools.

I. INTRODUCTION

The single crystal diffractometer (SCD) at the intense pulsed neutron source (IPNS) was the among first of its kind when it was built. It uses a four circle Huber diffractometer with the χ angle fixed at 45 degrees and a single angle camera at 90 degrees from the direct beam at 30 cm from the sample. The detector has 85 x 85 pixels in an active area of 30 cm x 30 cm.¹ This has been the instrument configuration, with minor adjustments, ever since. SCD is currently going through a major upgrade, replacing its original detector with two higher resolution angle cameras that each have an active area of 15 cm x 15 cm at a distance of 20 cm from the sample. In addition, the data acquisition system (DAS) for the SCD is also being upgraded. Not surprisingly the entire instrument upgrade is making it necessary to update much of the analysis software.

The Integrated Spectral Analysis Workbench (ISAW) is a general purpose data analysis and visualization package which can read the raw data of any IPNS instrument, including the upgraded SCD. ISAW also has an extendible system for performing calculations and operations on the data well suited to developing new analysis methods. This paper will describe the new data analysis package. Its early versions use much of the old software while implementing new methods of data visualization.

II. CURRENT SOFTWARE

The intention of all measurements on SCD is to collect up to 2 steradians of diffraction data from a single crystal sample. The raw data is then reduced to Bragg intensities consisting of the (hkl) index, integrated intensity, and uncertainty in the integrated intensity. This information is then used to determine the contents of the unit cell. To measure the full 2 steradians up to 44 crystal orientations need to be measured. For many crystals this is unnecessary and far fewer orientations are measured to reveal all of the unique reflections. An example of this

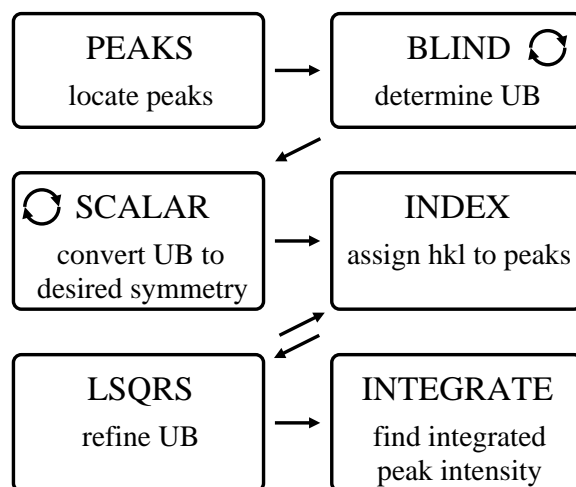


FIG. 1: Outline of the data analysis to produce a list of indexed, integrated peak intensities.

is quartz, with a hexagonal crystal structure, which only needs 13 crystal orientations measured to collect all of the unique reflections.

The original analysis software was written using FORTRAN 77 to run in a batch mode on the VAX. The analysis steps are outlined in Fig. 1. After the analysis is completed the list of reflections, frequently called peaks, are then converted to a format appropriate for the modeling software used. PEAKS determines the location of the brightest points in time and pixel space which is converted to position and wavelength. This is normally run once for each crystal orientation and all reflections are collected in the same file. Then BLIND is run repetitively with between 10 and 17 reflections to determine lattice parameters ($a, b, c, \alpha, \beta, \gamma$) and an orientation matrix (UB) which are satisfactory.^{2,3} SCALAR transforms the symmetry of the UB matrix found with BLIND to be more appropriate. For the quartz example mentioned above the UB matrix determined by BLIND will frequently make b the long axis, SCALAR performs

¹ A .Schultz et al., Transactions ACA 18, 169 (1982).

³ W .R .Busing and H .A .Levy, Acta Cryst. 22, 457 (1967).

² C .G iacovazzo et al., Fundamentals of Crystallography, Oxford Science Publications, 1992.

This figure "contour.jpg" is available in "jpg" format from:

<http://arxiv.org/ps/cond-mat/0210427v1>